


Sefar Solutions for Infusion Sets

SEFAR MEDIFAB fabric having a mesh aperture of < 20 µm ensures that during an infusion no foreign particles enter the patient's bloodstream. We offer tapes, tubes and bags.

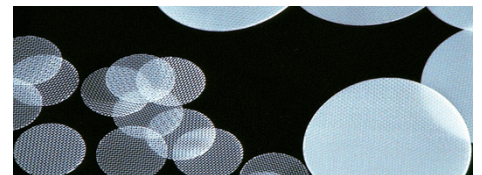

Product Features

Function

The infusion filters serve as safety filters. The function assures the retention of particles (ISO 8536-4). The source of such particles can be the infusion bottle, infusion bag or infusion solution (for example, aggregated nutrients or pharmaceutical compounds).

SEFAR MEDITEX membranes are used to ventilate infusion systems (ISO 8536-4). The construction of the membrane guarantees a high air permeability while preventing the ingress of viruses and bacteria as well as the escape of liquids.

DOWNLOADS

Brochure: Medical Components (PDF 3997 kb)

Customer Info: SEFAR MEDIFAB Fabrics are Low Protein Binding (PDF 208 kb)

Your Benefits

Fabrics

- Precise fabric dimensions
- Homogeneous material properties
- Exceptional lot-to-lot consistency
- Well-defined surface characteristics
- High throughput rates
- Excellent particle retention efficiencies

Sefar Inc.

111 Calumet Street
 Buffalo, NY 14043 – USA

Phone 800 995 0531
 Fax +1 716 685 9469

sales@sefar.us


[Go to product page](#)

- Non-shedding due to monofilament yarn
- Unsurpassed, two-dimensional stability
- Good sterilization


Membranes

- Homogenous material properties
- Exceptional lot-to-lot consistency
- High air permeability
- Minimal pressure loss
- Protection against viruses and bacteria
- Easy processing

Product range

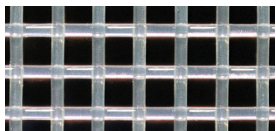

■ SEFAR MEDIFAB

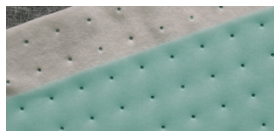

■ SEFAR MEDITEX
PORE


Sefar Inc.

111 Calumet Street
Buffalo, NY 14043 – USA

Phone 800 995 0531
Fax +1 716 685 9469

sales@sefar.us

[Go to product page](#)

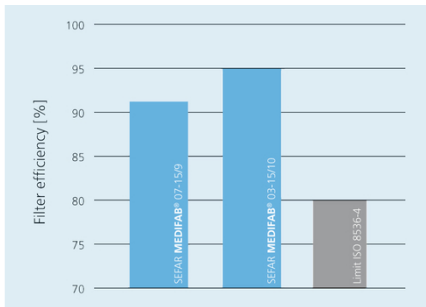

Filtration efficiency

Filtration efficiency of 20 µm particles of SEFAR MEDIFAB 03-15/10 and 07-15/9 in comparison to ISO requirements.

For further information and fabric samples, please select the preferred contact option on the right.

Locations



Sefar Inc.

Gray, GA – USA
 Phone: 800 995 0531
 Fax: +1 478 986 6953

[E-Mail](#)



Sefar Inc.

Kansas City, MO – USA
 Phone: 800 995 0531
 Fax: +1 816 452 2183

[E-Mail](#)



Sefar Inc.

Los Angeles, CA – USA
 Phone: 800 995 0531
 Fax: +1 909 544 5901

[E-Mail](#)



Sefar BDH Inc.

200 rue Clement Gilbert
Phone: +1 418 690 0888
Fax: +1 418 690 9499

■ [E-Mail](#)



Sefar Inc.

111 Calumet Street
Buffalo, NY 14043 – USA
Phone: 800 995 0531
Fax: +1 716 685 9469

■ [E-Mail](#)

Sefar Inc.

111 Calumet Street
Buffalo, NY 14043 – USA

Phone 800 995 0531
Fax +1 716 685 9469

sales@sefar.us

[Go to product page](#)